



**Thermo Scientific VersaCool**  
Refrigerated Circulating Bath



experience the benefits  
of going headless

**Thermo**  
SCIENTIFIC



# We lost our head, so you don't lose your cool


We didn't set out to design one of the world's most innovative refrigerated circulating baths. We simply set out to solve the issues that labs like yours face everyday – smaller workspaces, tighter budgets, laboratory stock issues, green-building goals, and overly complex equipment.

## Off with its head!

We started with the form factor: how to gain a large work area with a limited footprint. As with most great innovations – the solution was simple – we removed the control head to create more space. And we didn't stop there. Designed from the ground up – we added enhanced remote monitoring capabilities, energy efficiency, global voltage input, and other tools and features to give you control and optimal precision.

Imagine more space for your samples and smarter tools and features to simplify your daily tasks; so well designed anyone could use it.

Introducing the new Thermo Scientific™ VersaCool™ Refrigerated Circulating Bath. It's designed to help address your daily challenges, so you can stay focused on what you enjoy doing most – uncovering the answers to science's most perplexing questions.

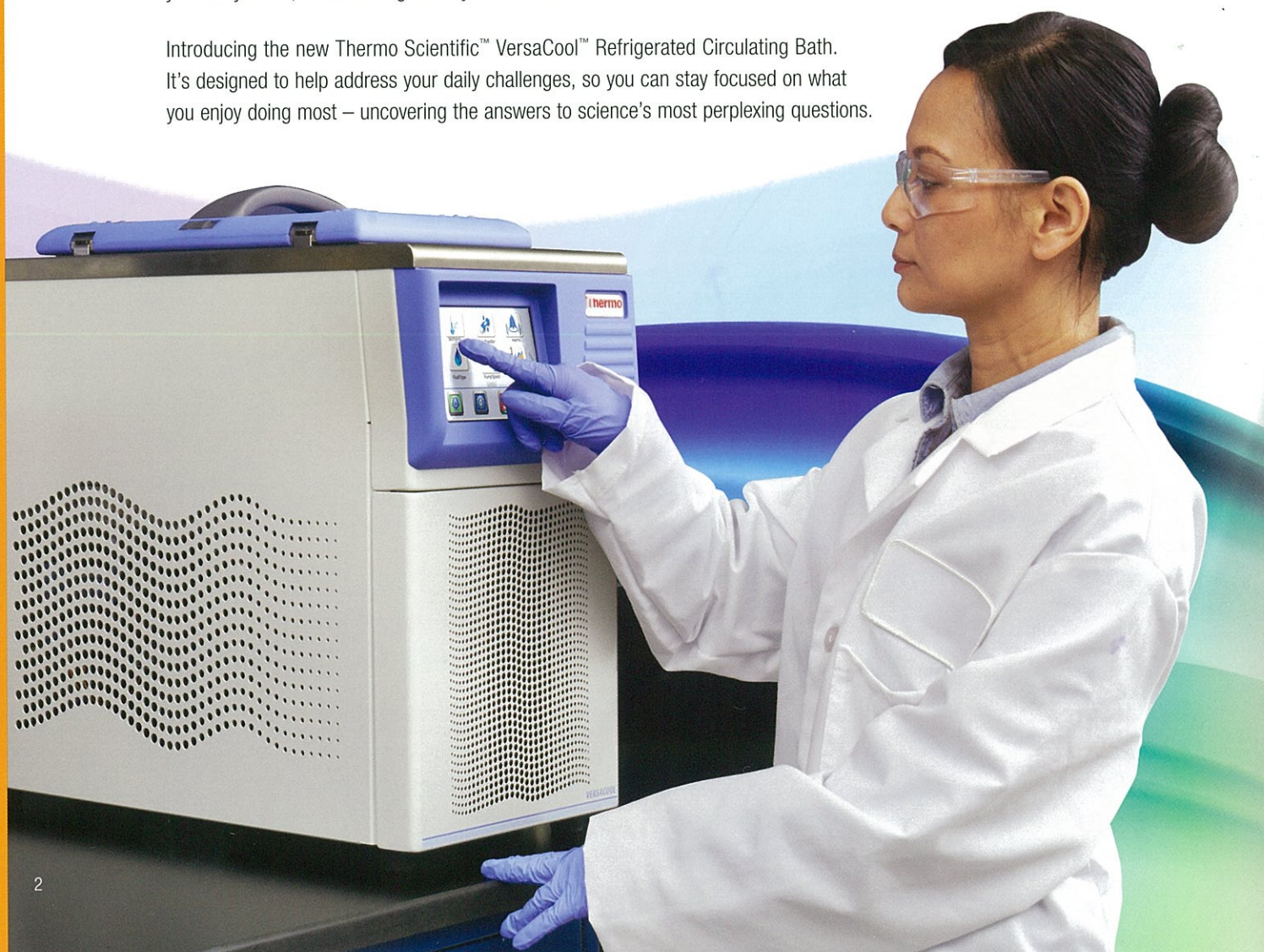





Redefining versatility,  
so you can  
accomplish more.

What does it take to be a leader?  
We believe it takes courage, confidence,  
insight, persistence, and perhaps the  
most important skill of all, versatility.

At Thermo Fisher Scientific, we understand  
the daily challenges and industry pressures  
that you may face. That's why we are  
continually evolving, innovating, and  
advancing our technologies to accelerate  
the path of scientific discovery and enhance  
value to customers.

That's our mission and promise to you.

VERSACOOOL




# Go headless, get more benefits



**We see our largest bath with the smallest footprint, you'll see more samples – one step closer to discovery**

- Removal of the control head and coils creates a larger, safer and easier to clean working area
- Increase capacity and the number of samples
- Maintain compact footprint on your workbench

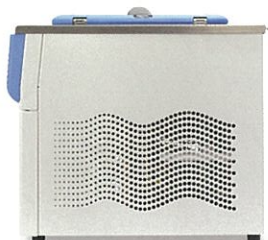

**We see our most energy efficient refrigerated circulating bath, you'll see lower power consumption and energy costs**

- Continuous operation in an energy efficient mode
- Variable speed control in pumps, fans and compressors
- Use only the energy you need for your application



**We see an intuitive, 5.7" touchscreen color display, you'll see an easier, more efficient set up**

- Glove/stylus friendly touchscreen interface
- Bright display to view critical readings
- Simple navigation and operation



**We see a tool-less, drip-less rack and lid, you'll see enhanced safety and faster cleanup**

- Tool-less rack accommodates a wide range of beaker and test tube sizes
- One lid that can be left-handed hinge, right-handed hinge, or unhinged
- Drip-less lid and optional gable cover reduce moisture and directs condensate back into bath



**We see global voltage input, you'll see one bath that works worldwide**

- Automatically detects the appropriate frequency and voltage
- One part number that works on all voltages worldwide
- Optimizes inventory management



**We see USB and Bluetooth™ communications, you'll see remote operation at your fingertips¹**

- Utilize the communication and control features to ensure your samples are properly maintained
- Pair to any iOS or Android smart phone or tablet for remote operation and monitoring
- Obtain operational control from remote locations with our optional NesCom software or with Thermo Scientific™ Smart-View™ wireless monitoring solution




# Think outside the box

Each VersaCool bath comes with external circulation connections so you can easily circulate the fluid from the bath to your application.




USB, MicroUSB, RS232 accessible from front of unit.<sup>2</sup>

<sup>2</sup> See specifications chart for details.

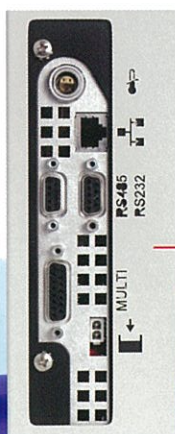

Comprehensive communications included.<sup>2</sup>

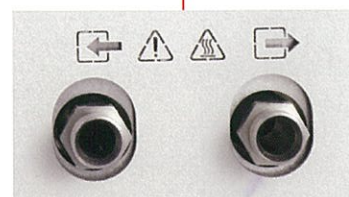


Hose barbs and threaded adapters included.<sup>2</sup>




## Performance curves

### Cooling Capacity

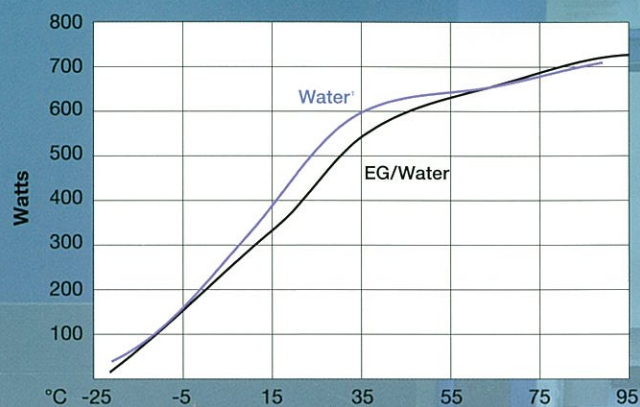

### Pump Performance

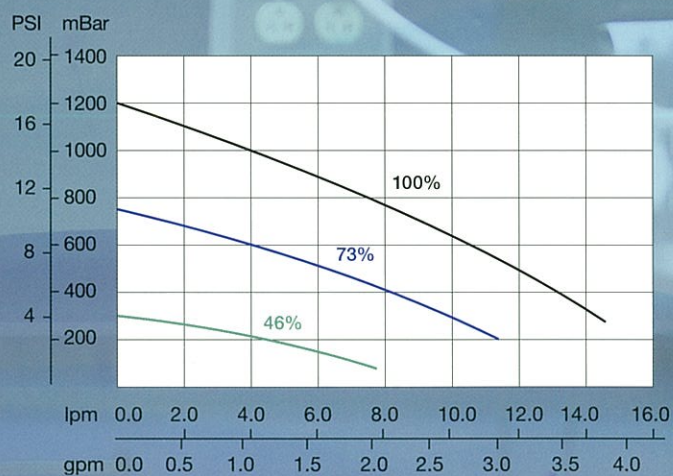

1. Specifications obtained at sea level using water (above 5°C to 90°C) or a fluid with a specific heat of 2.3 kJ/kg-K or 0.55 Btu/lb-F (less than 5°C) as the recirculating fluid at a 20°C ambient condition, at nominal operating voltage. Other fluids, process temperatures, ambient temperatures, altitude or operating voltage will affect performance. Specifications are for reference only and are subject to change.



# Frequently asked questions

## What's the difference between the VersaCool Refrigerated Circulating Bath and standard refrigerated circulating baths?

VersaCool Refrigerated Circulating Baths were designed with a focus on creating one of the largest work areas with the smallest footprint. Without the control head, heating element, cooling coils and pump in the bath area, you have more room to fit your samples, beakers, or test tubes. This headless design also makes it easier to clean and safer to operate. VersaCool also features a lid that can be configured to the user's comfort – left-handed hinge, right-handed hinge, or unhinged.

## How does Global Voltage work?

The VersaCool bath features Global Voltage for ordering and operational simplicity. It automatically detects the input voltage range with no user intervention for plug-and-play operation anywhere in the world. No matter where you or your customers are – one part number is all that is required.

## Does the VersaCool bath come with external circulation connections?

Yes! Not only is VersaCool equipped with one of the most powerful bath pumps in its class (with industry standard 16M x1 male connections), a complete selection of hose barbs and threaded adapters are included as standard.

- ¼", ½", 8mm and 12mm hose barbs
- ¼" MNPT and ½" MNPT – ideal for use with quick disconnects

## Does the VersaCool bath come with remote monitoring capabilities?

Yes, the VersaCool comes standard with RS232, RS485, USB, and MicroUSB for serial communication with a PC or application that can use those standards. Through these methods complete control, monitoring and data logging are available. VersaCool also has an Analog multifunction port that has contacts for fault/warnings detection, remote on/off and 4-20mA, 1mV/°C and 0-10VDC ranges for setting and reporting temperature to your PLC. Lastly VersaCool has Bluetooth™ and can be controlled and monitored using your Bluetooth enabled iOS or Android phone or tablet. An integrated probe well is located on the back of the unit to accommodate RTD monitoring probes including the optional Smart-View system probe.

## What programming parameters does VersaCool have?

VersaCool has a wide variety of operation, convenience and safety parameters that can be selected or programmed by the user. For temperature ramping profiles it includes the ability to program temperature profiles over time, chart and data log. You can save profiles and use them on multiple baths by utilizing the USB flash-drive download and upload capabilities. You can set temperature limits, alarms, pump speeds, and much more using the touchscreen interface.



# VersaCool Refrigerated Circulating Bath

Specifications	
Minimum temperature (°C)	-20
Maximum temperature (°C)	150
Cooling capacity (W)	425
Heater capacity (kW)	2 / 1.2
230V / 115V	
Maximum pressure (mbar / psi)	1200 / 17.4
Maximum flow rate (LPM / GPM)	14.6 / 3.86
Ambient temperature range (°C)	10 to 40
Noise (dBA)	<58
Voltage range (VAC / Hz)	100-120 / 50 or 60 200-240 / 50 or 60

## Dimensions / Weight

Working area dimensions D x W x L	15.0 x 17.3 x 29.7 (cm) 5.9 x 6.8 x 11.7 (in)
Overall dimensions H x W x L	57.0 x 25.8 x 47.4 (cm) 22.4 x 10.2 x 18.7 (in)
Weight (kg / lbs)	36 / 80

## Connectivity<sup>1</sup>

Bluetooth	✓
Ethernet	Pending
Micro USB – front	✓
USB – front	✓
RS232 – front / back	✓
RS485	✓
Remote sensor port	✓
Smart-View monitor port	✓

## Compliance

CE	✓
RoHS	✓
UL / CSA	Pending
WEEE	✓

## Warranty

3 / 2 warranty <sup>2</sup>	3 year warranty on electronic components 2 year warranty on all other components
-----------------------------	---

1. See product literature for technology requirements

2. Subject to the terms and conditions of the standard limited warranty found in Thermo Fisher Scientific's Standard Terms and Conditions of Sale, available at [www.thermoscientific.com](http://www.thermoscientific.com).

Ordering Information			
Part Number	Description	Part Number	Remote Temperature Sensors
231121300	VersaCool Refrigerated Circulating Bath (includes VersaLid & Plumbing Adapter Kit)	3330818	Pt100 probe, teflon coated, flexible, 300mm long, 3mm Ø, cable length 3m
230000001	Stainless Steel, Adjustable VersaRack • Holds up to 75 test tubes – 10mm Ø • Holds up to 48 test tubes – 16mm Ø • Holds up to 21 test tubes – 25mm Ø	3330429	Pt100 probe, 18/8 stainless steel tubing, 150mm long, 3mm Ø, cable length 3m, up to 600°C
Part Number	Heat Transfer Fluids		
610000000000	Silicone oil, temperature range +30°C to +150°C, 5 GAL		
610000000005	Algaecide / corrosion inhibitor, Nalco Kit		
610000000007	THERMO200 Treated water solution with Nalco, temperature range +5°C to +95°C, 5 GAL		
9990203	Sil 180 Silicone oil bath liquid, temperature range -40°C to +200°C, 5L		
9990204	Sil 180 Silicone oil bath liquid, temperature range -40°C to +200°C, 10L		
9990213	Synth 260 bath liquid, temperature range +40°C to +250°C, 5L		
9990214	Synth 260 bath liquid, temperature range +40°C to +250°C, 10L		
160000000001	Ethylglycol, 5 GAL (approx. 19L) for low temperature applications to -30°C		
Part Number	Software		
422000000004	NEScom control / monitoring PC software		
Part Number	Tubing and Accessories		
1600029	Adapts the M16x1 end of insulated stainless steel flexible tubing to 1/4" male pipe thread (MNPT). Contains 2 adapters		
1600146	Plumbing Package – includes (4) clamps and (2) 5' Viton tubing (uninsulated), temperature range of -30°C to +200°C, 12mm Ø		
1600147	Plumbing Package – includes (4) clamps and (2) 5' Viton tubing (insulated), temperature range of -30°C to +200°C, 12mm Ø		
3330293	Insulated flexible stainless steel tubing with M16x1 threaded fittings 1.0m long. -50°C to +300°C		
3330292	Insulated flexible stainless steel tubing with M16x1 threaded fittings 0.5m long. -50°C to +300°C		
3330294	Insulated flexible stainless steel tubing with M16x1 threaded fittings 1.5m long. -50°C to +300°C		
0012560	Connects two of the insulated flexible stainless steel tubing together		

## Add Thermo Scientific Smart-View wireless remote monitoring<sup>3</sup>

Safeguard the integrity of your precious samples by monitoring temperature and securely logging data to give you peace of mind.

Visit [www.thermoscientific.com/smart-view](http://www.thermoscientific.com/smart-view) to learn more.

3. Smart-View is not available in all Radio Frequency (RF) regions. Please ask your sales representative to evaluate a Smart-View wireless solution that best meets your specific needs.



Experience the benefits of going headless at [thermoscientific.com/versacool](http://thermoscientific.com/versacool)

© 2014 Thermo Fisher Scientific Inc. All rights reserved. Bluetooth is a registered trademark of SIG, Android is a registered trademark of Google. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Results may vary under different operating conditions. Specifications, terms and pricing are subject to change. Smart-View is not available in all countries. Please consult your local sales representative for details.

**North America:** USA/Canada tollfree: +1 (800) 258-0830; USA: +1 (603) 436-9444 or [info.tc.us@thermofisher.com](mailto:info.tc.us@thermofisher.com)

**Europe:** Benelux: +31 (0) 76 579 55 55 or [info.tc.nl@thermofisher.com](mailto:info.tc.nl@thermofisher.com); France: +33 (0) 1 60 92 48 00 or [info.tc.fr@thermofisher.com](mailto:info.tc.fr@thermofisher.com);

Germany: +49 (0) 721 4 09 44 44 or [support.tc.de@thermofisher.com](mailto:support.tc.de@thermofisher.com); United Kingdom: +44 (0) 8706 099 254 or [info.tc.uk@thermofisher.com](mailto:info.tc.uk@thermofisher.com)

**Asia:** China: +86 (21) 68 65 45 88 or [info.tc.china@thermofisher.com](mailto:info.tc.china@thermofisher.com); India: +91 (22) 27 78 11 01 or [info.tc.in@thermofisher.com](mailto:info.tc.in@thermofisher.com);

Japan: +81 45 453 9220 or [info.lpg.jp@thermofisher.com](mailto:info.lpg.jp@thermofisher.com)

BRTCVERSACOOOL 0314

**Thermo**  
SCIENTIFIC

Part of Thermo Fisher Scientific