


## High Performance UV Transilluminators

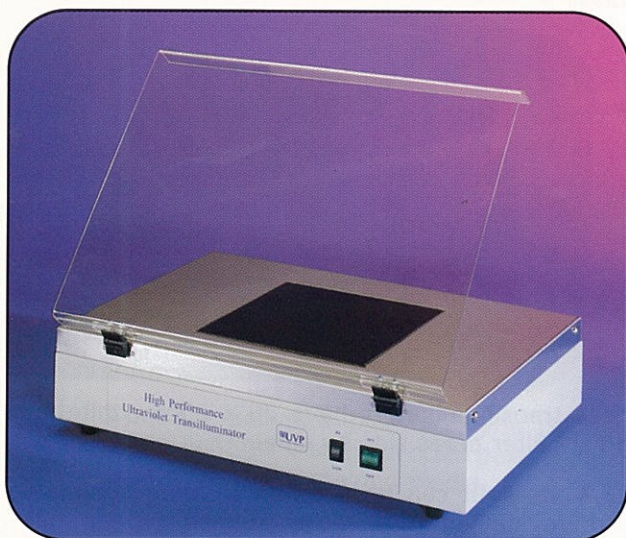

**High Performance Transilluminators**  
Hinged cover is adjustable to different angles.

The exclusive 25-watt ultraviolet tubes provide a total of 100-watts of high intensity performance and brilliant ultraviolet illumination.

- Select from shortwave (254nm), midrange (302nm) or longwave (365nm) UV for high fluorescence and sensitivity on stained gels
- Stainless steel frame enables easy cleaning
- Brilliant back-lit UV fluorescent illumination is enhanced with UVP's UVG long-life filter
- The transilluminator's powerful cooling fan maintains a low filter surface temperature to prevent sample damage and maintain a constant UV intensity
- Removable and adjustable protective UV blocking cover is included with each transilluminator
- Models with Hi/Lo intensity settings enable the researcher to switch settings depending on research requirements

Dimensions:

13.25D x 19.13W x 5.63H in. (337 x 486 x 143mm)

Height includes UV blocking cover.

Weight: 21 lb. (9.6 kg)

### Ordering Information

<b>TFM-20 Hi/Lo Intensity</b> 302nm UV Filter Size 20 x 20cm Replacement Filter	95-0286-01 (115V) 95-0286-02 (230V)  38-0190-02
<b>TFM-26 Hi/Lo Intensity</b> 302nm UV Filter Size 21 x 26cm Replacement Filter	95-0285-01 (115V) 95-0285-02 (230V)  38-0190-09
<b>TFM-30 Hi/Lo Intensity</b> 302nm UV Filter Size 25 x 30cm Replacement Filter	95-0302-01 (115V) 95-0302-02 (230V)  38-0190-11
<b>TFM-40 Hi/Lo Intensity</b> 302nm UV Filter Size 20 x 40cm Replacement Filter	95-0283-04 (115V) 95-0283-05 (230V)  38-0190-04
<b>TFL-40 Hi/Lo Intensity</b> 365nm UV Filter Size 20 x 40cm Replacement Filter	95-0283-01 (115V) 95-0283-02 (230V)  38-0190-10
<b>TFS-20 Hi/Lo Intensity</b> 254nm UV Filter Size 20 x 20cm Replacement Filter	95-0319-01 (115V) 95-0319-02 (230V)  38-0190-02
<b>TFS-26 Hi/Lo Intensity</b> 254nm UV Filter Size 21 x 26cm Replacement Filter	95-0320-01 (115V) 95-0320-02 (230V)  38-0190-09
<b>TFS-30 Hi/Lo Intensity</b> 254nm UV Filter Size 25 x 30cm Replacement Filter	95-0321-01 (115V) 95-0321-02 (230V)  38-0190-11
<b>TFS-40 Hi/Lo Intensity</b> 254nm UV Filter Size 20 x 40cm Replacement Filter	95-0322-01 (115V) 95-0322-02 (230V)  38-0190-04

### Replacement Tubes

Tube, 25-watt, 254nm	34-0073-01
Tube, 25-watt, 302nm	34-0072-01
Tube, 25-watt, 365nm	34-0060-01

## 2UV High Performance Transilluminators

2UV models are the same size and design as other high performance models but are equipped with midrange (302nm) and longwave (365nm) UV in each unit.

### Ordering Information

<b>TFML-26E Single Intensity</b> 302/365nm UV Filter Size 21 x 26cm Replacement Filter	95-0308-01 (115V) 95-0308-02 (230V)  38-0190-09
<b>TFML-40E Single Intensity</b> 302/365nm UV Filter Size 20 x 40cm Replacement Filter	95-0309-01 (115V) 95-0309-02 (230V)  38-0190-04
<b>TFML-20E Single Intensity</b> 302/365nm UV Filter Size 20 x 20cm Replacement Filter	95-0323-01 (115V) 95-0323-02 (230V)  38-0190-02
<b>TFML-30E Single Intensity</b> 302/365nm UV Filter Size 25 x 30cm Replacement Filter	95-0324-01 (115V) 95-0324-02 (230V)  38-0190-12