


Grinding/Polishing



MODEL: P20F-1 P25F-1

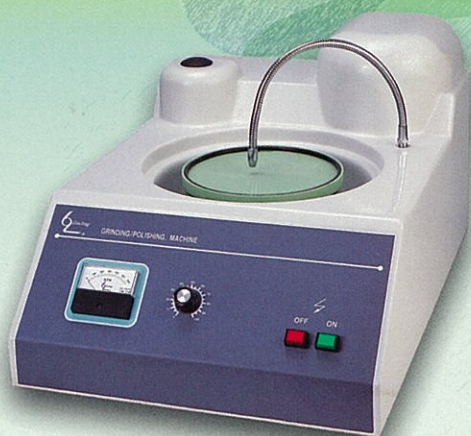

MODEL: P20FR P25FR

Features

1. Precise rotational balance will not deviate even heavily pressed.
2. Steady machine design does not waver during operation.
3. It is designed with double waterproof lids, safe and easy grinding, which makes the workpiece uneasy to fall.
4. The main shaft is anti-leakage designed, which makes running smooth.
5. The replaceable grinding disc can be changed quickly.
6. Double stainless steel water pipes can be adjusted at will.
7. Special O ring rubber design fixed sand paper and polishing flannel prevent press lids from ejection and hurting hands.
8. Workpiece rotates regularly with fast and neat grinding.

MODEL	P20F-1	P25F-1	P20FR	P25FR
No. of discs	1		2	
Grinding Method	Manual		Manual	
Disc dia	φ 200	φ 250	φ 200	φ 250
Rotational speed	0 ~ 600rpm		0 ~ 600rpm	
Motor power	DC180V or DC90V 1φ 180w		DC180V or DC90V 1φ 250w	
Cooling unit	1 piece		2 piece	
Machine type	Table top		Table top	
Dimensions	740(W) x 440(D) x 360(H)mm		740(W) x 700(D) x 360(H)mm	
Power	AC 220V 1φ or AC 110V 1φ		AC 220V 1φ or AC 110V 1φ	
Weight	28kg	30kg	45kg	50kg
Standard accessories	Aluminum Oxide Powder 100g		Aluminum Oxide Powder 100g	
	Drain Pipe 1.5M	Anti-splash Cover 1pc	Drain Pipe 1.5M	Anti-splash Cover 2pc
	Small Pipe Clip 2pcs	Water Sand Paper 1pc	Small Pipe Clip 2pcs	Water Sand Paper 2pc
	Abrasive Cloth 1pc	Big Pipe Clip 2pcs	Abrasive Cloth 2pc	Big Pipe Clip 2pcs
	Water Inlet Pipe 1.5M	O Ring 1pc	Water Inlet Pipe 1.5M	O Ring 1pc